

# SmileGoods®


## CORDLESS WATER FLOSSER

### QUICK GUIDE AND USER MANUAL



## SYSTEM COMPONENTS

Item	Description	Qty
A	Water Flosser Power Handle	1
B	Standard Jet Tip	2
C	Periodontal Pocket Tip	1
D	Orthodontic Tip	1
E	USB Charging Cable	1
F	Quick Guide and User Manual	1

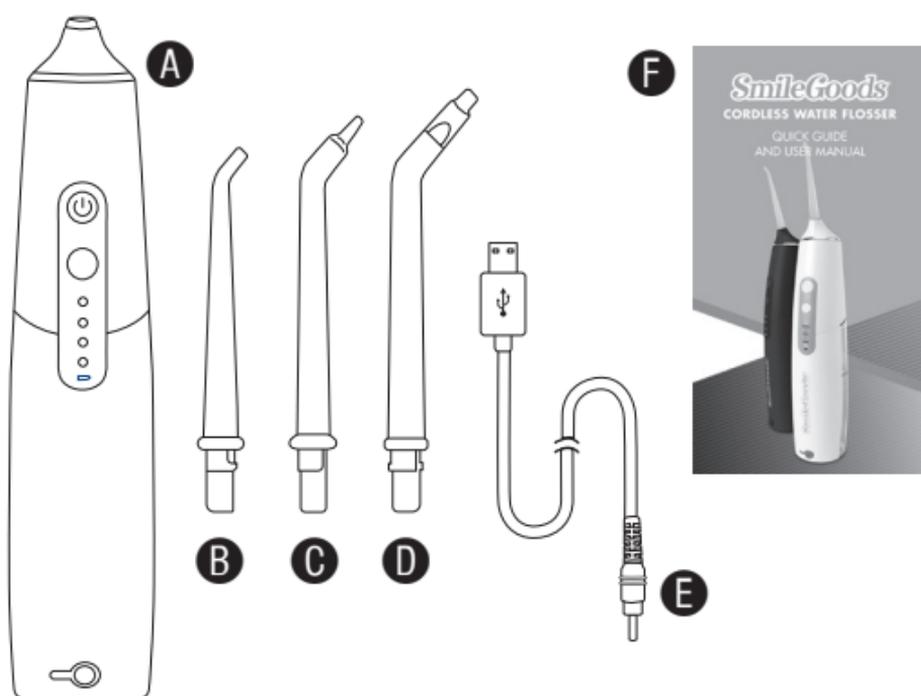

# PRACTICON

PRACTICAL INNOVATIONS FOR DENTISTRY  
1-800-959-9505 • Fax 1-800-278-0885  
practicon.com • Greenville, NC 27834  
Made in China • IFU 7155012\_1

# **THANK YOU** for purchasing the SmileGoods Cordless Water Flosser. With proper use and care, the Cordless Water Flosser will help you maintain a healthy smile for years to come.

This product is designed for cleaning between teeth, beneath prosthodontic restorations, and around orthodontic appliances. It is also intended for irrigating and flushing periodontal pockets and intraoral surgical sites to aid in healing and pain relief. First-time users may experience minor gum irritation or bleeding. If condition persists for more than 1–2 weeks, discontinue use and consult a dental professional.

Use only as directed. Prosthodontic, orthodontic, surgical patients, and patients under 6 years of age should use under the guidance of a dental professional.

In order to reduce the risk of failure, damage, fire, electric shock, injury or other hazards, please read and understand all instructions before using this device. Save this manual for your reference.

## **GETTING STARTED**

To use, select appropriate nozzle tip for your health needs. Choose standard jet tip for most conditions. For periodontal irrigation and sensitive gums, use the periodontal pocket tip. The orthodontic tip is specifically designed to clean around orthodontic wires and brackets.

Fully insert the tip into the power handle, and turn tip 40 degrees clockwise to lock into place. To remove tip, turn it counterclockwise and pull out gently. Do not use power handle without a tip installed.



To fill, open the reservoir door and hold handle horizontally while filling. Close door securely. Do not fill reservoir with mouthwash or a strong antibacterial agent without proper dilution.

Check for proper tip installation and insert tip into mouth **BEFORE** pressing the power button. When using, aim the nozzle perpendicular to the irrigation site and keep mouth slightly open.

**POWER BUTTON**

**MODE/PRESSURE BUTTON**

## **THREE OPERATING MODES**

**Hi/Lo NORMAL** » Highest pressure for deep cleaning and regular flossing

**Hi/Lo SOFT** » Gentle pressure for first time users and sensitive gums

**Hi/Lo PULSE** » Ideal for gum massage and intraoral irrigation

**CHARGING PORT**

## **THREE-MODE OPERATION**

*with Hi and Lo Power Settings*

The SmileGoods Cordless Water Flosser provides three operating modes with Hi and Lo pressure settings for all user preferences and oral health needs.

To turn ON, press the power button. The unit will operate for 2 minutes and turn OFF automatically. To stop the cleaning at any time during a cycle, press the power button again.

Press the Mode/Pressure button repeatedly to cycle through the six available user settings indicated by the corresponding LED position and brightness.

2min

**2-minute shut-off timer**

## **USB CHARGING**

The SmileGoods Cordless Water Flosser is rechargeable with USB cable. Before first use and at least once every 3 months thereafter, charge unit for a full 24 hours.

*To charge:*

- Connect the charging cable to a USB power source (5V@1A minimum, not included).
- Insert cable into charging port at base of power handle.
- While charging, battery LED indicator will change color to show charging status:
  - Red <20% charge
  - Blue 20% to 100% charge
  - Green 100% charge

During use, battery LED indicator will change color to show battery status:

- Green 100% to 90% charge
- Blue 90% to 20% charge
- Red <20% charge

After each use, it is safe to store your flosser with the USB charger connected, if preferred. Once per year, allow battery to drain completely before recharging.



## CARE AND CLEANING

### Nozzle tips

Rinse with water and wipe as needed. Do not bend or alter the tips. For best results, replace tips every 6 months or if damaged.

### Power handle

Rinse with cold water as needed and wipe dry. Do not use bleach or acid-based cleaners. Do not immerse the power handle in water or any other liquid. If the unit will not be used for 1 month or more, empty the reservoir and store with fill door open.

## Q&A

### Q: Why is the water pressure getting weaker?

A: The battery may be low and require charging, or the tip could be blocked or installed incorrectly. Be sure tip is inserted correctly. Check for debris inside tip and replace tip if damaged.

### Q: Why is the battery life still short after fully charging?

A: Ensure water flosser charges for a full 24 hours before next use, and that the charging cable is fully inserted into charging port. Check that USB power source is sufficient (5V@1A minimum).

### Q: After long periods without use, why is the battery indicator light not working?

A: When not used for 3 months, the battery will enter auto-protection mode and require a 24-hour charge before use again.

## SPECIFICATIONS

Item: 7155012

Weight: 11.2 oz.

Material: ABS and Polycarbonate

Modes: Hi/Lo @ Normal/Soft/Pulse

Water pressure: 30–110 psi

Reservoir capacity: 6.8 fl. oz. (200ml)

Battery: 2,000mAh Lithium-ion battery

Charge time: 4 hrs.

Battery life: 60 mins.

Waterproof rating: IPX7 (1m @ 30 mins)

Certificate: FCC, FDA, RoHS, IPX7

Warranty: 1 year



## WARNINGS

### Electric Shock Hazard

The **SmileGoods Cordless Water Flosser** is an electric device designed to meet worldwide safety standards.

To avoid electric shock:

- Do not attempt to open or alter the unit in any way. Any attempt to open the unit will void the warranty.
- Do not put any foreign or metal objects into the reservoir or the unit housing.

### Disposal

Please dispose of unit according to local environmental guidelines or regulations regarding lithium-ion batteries. DO NOT dispose of unit in fire.

### Battery

Do not charge, operate, or store at temperatures above 104°F (40°C). If unit is stored for 3 months or longer, recharge battery for 24 hours before using.

### Repairs and Defects

Do not use if the unit is damaged or defective. If battery is defective, do not transport the unit by air under any circumstances.

Only use **SmileGoods Cordless Water Flosser Replacement Tips** available by calling 1-800-959-9505 or at [practicon.com](http://practicon.com).

Item Number: 71-55013

Description: 4 replacement tips (2 standard jet, 1 periodontal pocket, 1 orthodontic)

## WARRANTY

Practicon, Inc. warrants this product to be free of manufacturing defects for a period of one (1) year from the date of purchase, deemed as the date of the invoice. At Practicon's discretion, units returned due to manufacturing defects within the warranty period will be repaired or replaced. The following are expressly excluded from the warranty:

- Damage and/or failure of the unit caused by dropping or mishandling.
- Damage and/or failure of the unit caused by natural disasters, such as earthquakes, floods, lightning, and voltage spikes.

Any attempts to open or repair the unit will invalidate the warranty.

Date of invoice \_\_\_\_\_